




Ref. 5021- 5022

GENERAL CHARACTERISTICS


Handling simplicity is what defines this line of autoclaves in its two versions: 11 litres and 22 litres. They comply with the standards of instrument and appliance manufacturers.


- Five programs: 4 sterilisation programs and 1 disinfection program.
- LCD screen controlled by a microprocessor, which indicates the real time of any process as well as any incidents that might occur.
- Program readout in six languages.
- Dual distilled water tank that prevents the recirculation of distilled water.
- Pre-vacuum and post-vacuum pump that prevents air bubbles for complete sterilisation and a subsequent vacuum for complete drying, free from bacteria using a bacteriological filter.
- Cycle duration:
 - 134° C: **3 minutes** -134° C: 15 minutes / 12 minutes of drying - **121° C: 30 minutes / 12 minutes of drying**
 - 121° C: 13 minutes** -105° C: 30 minutes / disinfection.
- Trays
- Tray holder
- Cassette tray
- Drain hose
- Tray removal handle
- Adapter for PC printer.


Manufactured according to EC directives

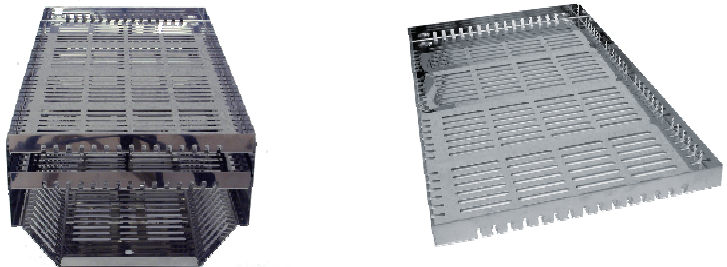

Ref.	Capacity	Dimensions W x H x D	Weight.	Consumption	Chamber	Chamber depth
5021	11 L.	480 x 380 x 520mm.	38 Kg.	660 W.	Ø200 mm.	340 mm.
5022	22 L.	530 x 430 x 600mm.	56 Kg.	900 W.	Ø250 mm.	470 mm.

Basket	Sheet blind basket		
		Ø x Long	Ref.
	Inox	Ø220x220mm.	CE 001
	Inox	Ø280x280mm.	CE 002
	Inox	Ø380x380mm.	CE 003
	Inox	Ø280x200mm	CE 004

Basket	Perforated basket		
		Ø x Long	Ref.
	Inox	Ø220x220mm.	CE 005
	Inox	Ø280x280mm.	CE 006
	Inox	Ø380x380mm.	CE 007
	Inox	Ø280x200mm	CE 008

Basket	Accessories information		
		Ø x Long	Ref.
	Inox	Ø220x220mm.	CE 009
	Inox	Ø280x280mm.	CE 010
	Inox	Ø380x380mm.	CE 011
	Inox	Ø280x200mm	CE 012