



www.bunsen.es



TUBE MIXERS

Ref. 1100/1110

GENERAL CHARACTERISTICS

- * Head stock for mixing test tubes of up to Ø 20 mm.
 - * It has two switchable working modes:
 - Continuous mode for uninterrupted mixing.
 - Manual mode in which mixing is activated automatically by slight pressure from the tube on the head stock.
 - * Metal housing painted in epoxy.
 - * Electronically regulated speed control.

Manufactured according to CE directives.




BUNSEN SA. C/ Río Alberche, 30, - 28970 - Humanes de Madrid, Madrid. SPAIN
Tlf. + 34 916 113 584 - Fax. + 34 916 128 254
info@bunsen.es - www.bunsen.es



Ref	Mixing speed (rpm)	Power (W)	Dimensions (W x H x D, mm)	Weight (Kg)	Eccentricity
1100	Variable between 50 and 2600	10	140 x 100 x 155	4	6 mm.
1110	Fixed 2600	10	140 x 100 x 155	4	6 mm.

Ref	Accessories
1110-1	Tube rack for 16 test tubes (Ø 15 mm max.). CONSULT
1110-2	Tube rack for 25 test tubes (Ø 12 mm max.). CONSULT
1110-3	Head stock for 6 Eppendorf tubes
1110-4	90-mm platform for containers of up to 250 ml

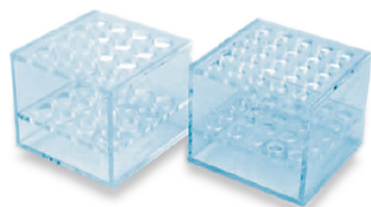

BUNSEN SA. C/ Río Alberche, 30, - 28970 - Humanes de Madrid, Madrid. SPAIN
Tlf. + 34 916 113 584 - Fax. + 34 916 128 254
info@bunsen.es - www.bunsen.es

