



www.bunsen.es

High SPEED



BLOBEL CENTRIFUGE

Ref. 1467

GENERAL CHARACTERISTICS

- * Maintenance – free, brush less motor.
- * Controlled by microprocessor with large LCD screen that permits displaying all parameters. 99 operation programs.
- * Timer: 0 a 99 hours 59 min and 59sec. or ∞ . (1 sec interval adjustment). Time counting begins either when push button start is pressed or when rpm of the rotor reaches rotational speed set up in the program.
- * Rotor identification. Automatic lid opening. Spinning in “short” mode.
- * Pre-selection for rotational speed of RCF. Possibility of varying parameters during the cycle.
- * Possibility of centrifuging samples with a density higher than 1.2gr / cm³ with automatic speed adjustment.
- * Centrifugal radius adjustment for samples with automatic RCF correction.
 - * Centrifuging parameters recording - PC computer (USB).
 - * Selected functions blocking, password protection.
- * Centrifugal radius adjustment for samples with automatic RCF correction.
 - * Menu language setting.



BUNSEN SA. C/ Rio Alberche, 30 - 28970 Humanes de Madrid, Madrid, SPAIN
Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es



- * Errors indication system: Signals of failure in the power supplier. Improperly chosen rotor.
- * Safety features:
 - * Lid locking during rotor running.
 - * Blocking at opened lid.
 - * Emergency lid lock release.
 - * Unbalance sensor.

Manufactured according to National and International standards EN-61010-1 and
EN-61010-2-020 on safety and EC marking.



Ref	Dimensions W x H x D	Weight	Power	Max. Capacity	Max. Speed	Max. RCF
1467	515 x 455 x 650mm.	74 Kg.	750W.	2000ml.	100 ÷ 18000 Rpm.	31150 (xg)
Timer (min)	Voltage		Display	Rotor ID	Refrigeration	Temperature Range
0-99h.-inf	115 v/ 230 v -- 50/60Hz		LCD	Si	No	-----

Ref	Accessories
1467-1761	 Angular 30 x 1.5/2.2 mL. / Max: 16400 Rpm. / RCF Max: 29198 xg 45°. Hermetically sealed lid.
1467-1762	 Angular 36 x 1.5/2.2 mL. / Max: 18000 Rpm. / RCF Max: 30065 xg 45°. Hermetically sealed lid.
1467-1763	 Angular 48 x 1.5/2.2 mL. / Max: 16400 Rpm. / RCF Max: 28265 xg 45°. Hermetically sealed lid.
1467-1766	 Angular 12 x 10 mL. / Max: 18000 Rpm. / RCF Max: 31150 xg / 30°. Hermetically sealed lid.
1467-1770	 Angular 12 x 15 mL. / Max: 14000 Rpm. / RCF Max: 23227 xg / 45°. Hermetically sealed lid. Falcon Tubes.
1467-1772	 Angular 8 x 30 mL. / Max: 17500 Rpm. / RCF Max: 29787 xg / 30°. Hermetically sealed lid. Nalgene type tubes.



BUNSEN SA. C/ Rio Alberche, 30 - 28970 Humanes de Madrid, Madrid, SPAIN
 Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es

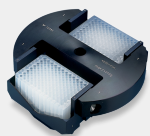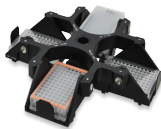

Ref	Accessories	
1467-1773		Angular 6 x 50 mL. / Max: 14500 Rpm. / RCF Max: 21625 xg / 30°. Hermetically sealed lid. Nalgene type tubes.
1467-1775		Angular 8 x 50 mL. / Max: 14000 Rpm. / RCF Max: 21255 xg / 30°. Hermetically sealed lid. Nalgene type tubes.
1467-1776		Angular 8 x 50 mL. / Max: 14000 Rpm. / RCF Max: 23666 xg / 30°. Hermetically sealed lid. Falcon Tubes.
1467-1777		Angular 4 x 250 mL. / Max: 10000 Rpm. / RCF Max: 14086 xg / 25°. Hermetically sealed lid.
1467-1778		Angular 6 x 100 mL. / Max: 12000 Rpm. / RCF Max: 17710 xg / 30°. Hermetically sealed lid. Falcon Tubes.
1467-1788C		Angular 36 x 10/15 mL. / Max: 5000 Rpm. / RCF Max: 3996 xg / 30°. Falcon Tubes.
1467-1789C		Angular 48 x 10/15 mL. / Max: 4700 Rpm. / RCF Max: 4025 xg / 30°. Falcon Tubes.
1467-2787C		Swing out 2 x 3 Microtiter or Blocks/ Max. 4300 rpm. / RCF Max. 2397 xg.
1467-2788C		Swing out 4 x 3 Microtiter or Blocks / Max. 4000 rpm. / RCF Max. 2826 xg.
1467-2786		Swing out 4 x 250 mL. / Max. 5000 rpm. / RCF Max. 4724 xg.
1467-2786C		Swing out 4 x 250 mL. / Max. 5000 rpm. / RCF Max. 4724 xg. Biológico safety.

Ref	Accesorios	
1467-3180		Swing out 8 x 50 ml. / Max. 5200 rpm. / RCF Max. 4413 xg. Falcon tubes.
1467-2830		Swing out 4 x 500 ml. / Max. 4300 rpm. / RCF Max. 3824 xg.
1467-2900		Swing out 4 x 750 ml. / Max. 4100 rpm. / RCF Max. 3458 xg.

Más centrifugas



Other centrifuges



BUNSEN SA. C/ Rio Alberche, 30 - 28970 Humanes de Madrid, Madrid, SPAIN
Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es

