



www.bunsen.es

DRYING OVEN



180L. DRYING OVEN WITH FORCED AIR CONVECTION Ref. 8328 Plus

GENERAL CHARACTERISTICS

- * Interior and exterior construction in stainless steel AISI 304.
- * Solid door, fitted with a key. Excellent thermal insulation.
- * Protection against overtemperature class 3.1 according to DIN 12880.
 - * USB port. Resolution 0.1°C.
- * PID microprocessor controller with a large (5,7") full colour touch screen, intuitive menu and user friendly software. They can be connected to Ethernet network for remote control.
- * Administrator function to manage User accounts
- * Access control via login. 7 days programming.
- * multi-segment temperature-time profile (up to 100). Loop function up to 99 times or endless.
- * Adjustable start delay feature (from 1 min to 99:59 h).
- * Overview of set and current parameters while operating. Adjustable ramps.
- * Recording of min, average and max temperature value for each segment.
- * Possibility of temperature calibration by the user.
- * Auto-diagnostic function. Automatic air-flap control.
- * Automatic fan shut-down after completing the program.
- * Audible and visual temperature alarm.
- * Power failure control system (program continued after restoring power).



BUNSEN SA. C/ Tokio, 2 28970 Humanes de Madrid, Madrid, SPAIN
Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es




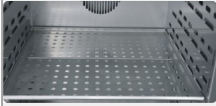






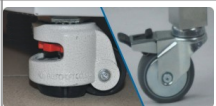
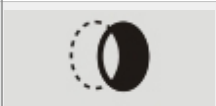
* Trays made of stainless steel wire; 3 tray included, maximum capacity; 9. Max shelf workload 25 kg. Consult accessories.

* Real time clock.

* Internal memory to store up to 4100 data records. And 20 user programs memory with password protected settings.

Manufactured according to CE directives.

| Ref. | Capacity | Dimensions W x H x D | Useful measures W x H x D | Weight | Temperature | Nominal power |
|--------------|----------|----------------------|---------------------------|--------|-------------|---------------|
| 8328P | 180L. | 650 x 1030 x 820mm. | 470 x 720 x 560mm. | 94 Kg. | Up to 300°C | 2800W |

| Ref | Accesories | |
|----------------|---|--|
| 8300-02 |  | Stainless steel wire shelf INOX |
| 8300-03 |  | Perforated shelf |
| 8300-04 |  | Reinforced shelf |
| 8300-05 |  | Reinforced version |
| 8300-06 |  | Interior led lighting |
| 8300-07 |  | Table with wheels |
| 8300-08 |  | Additional Pt 100 temperature sensor |
| 8300-09 |  | HEPA-fresh air filter |
| 8300-10 |  | Wheels |
| 8300-11 |  | Non-standard access port for external sensor |