



www.bunsen.es

Microcentrifuges



TATTUS S CENTRIFUGE

Ref. 1480s

GENERAL CHARACTERISTICS

Small laboratory centrifuge. Its structure guarantees work and the wide range of applications in medicinal, veterinary, biochemical and other laboratories where constituents of investigated material are to be separated under the effect of centrifugal force.

- * Maintenance – free brushless motor.
- * Short time of acceleration.
- * Controlled by microprocessor with LCD screen that permits displaying all parameters.
 - * Programmable in either RPM or RCF without distinction.
 - * Inspection glass.
- * Safety features:
 - * it is not possible to open the lid during spinning.
 - * it is not possible to start when the lid is open.
 - * Emergency lid lock release.

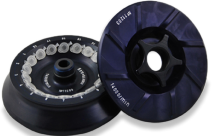
Manufactured according to National and International standards EN-61010-1 and EN-61010-2-020 on safety and EC marking.



BUNSEN SA. C/ Rio Alberche, 30 - 28970 Humanes de Madrid, Madrid, SPAIN
Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es



Ref	Dimensions W x H x D	Weight	Power	Max. Capacity	Max. Speed	Max. RCF
1480 S	220 x 180 x 270mm.	5 Kg.	75W.	90 mL.	100 / 6000 Rpm.	3340 (xg)
Timer (min)	Voltage		Display	Rotor ID	Refrigeration	Temperature Range
0-60 -inf	115 v/ 230 v -- 50/60Hz		LCD	No	No	-----

Ref	Accessories	
1480 1140C		Angular 6 x 10/15mL. / Max: 5800 Rpm. / RCF Max: 3122 xg / 30°
1480 1141C		Angular 8 x 10/15mL. / Max: 5800 Rpm. / RCF Max: 3122 xg / 30°
1480 1202		Angular 12 x 1.5/2.2mL. / Max: 6000 Rpm. / RCF Max: 2616 xg / 45° Hermetically sealed lid
1480 1204		Angular 24 x 1.5/2.2mL. / Max: 6000 Rpm. / RCF Max: 2616 xg / 45° Hermetically sealed lid
1480 1203		Angular 18 x 0.5mL. / Max: 6000 Rpm. / RCF Max: 2375 xg / 45° Hermetically sealed lid

Más centrifugas



Other centrifuges



BUNSEN SA. C/ Rio Alberche, 30 - 28970 Humanes de Madrid, Madrid, SPAIN
Tlf. +34 916 113 584 - Fax. +34 916 128 254
info@bunsen.es - www.bunsen.es

